



safety first 



DRIVE IN STYLE

The Permlock C3 allows you to drive or ride while sitting in your Permobil. This device secures your wheelchair inside the vehicle and is designed for use with the Permobil C300 and C350 wheelchair bases. You can adjust both the angle and height of the Permlock C3, which helps when mounting it in a vehicle with a sloped driver's position. A locking device mounted to the wheelchair chassis and a locking plate that is mounted to the vehicle ensure the wheelchair will not break free in the event of a collision or other sudden movement. In addition, you can pull either forwards or backwards into the locking plate.

**PERMOLOCK
C3**



SPRING-LOADED PINS

Two fixed pins are attached to the wheelchair chassis, which push up into the chassis when an obstacle is encountered, ensure a safe and secure ride. When an obstacle is encountered, they push up into the chassis.



LOCKING PLATE

Easy-to-install docking station that holds the chair firmly and securely in position.



ADJUSTABLE BASE

The height of the base can be adjusted in 5mm increments, and the angle adjustment utilizes a spring to help mount it on a slope.



VEHICLE CONTROL BUTTON

Located on the dash, it has three lights, -red, amber and green, to indicate if the chair is locked into place. Simply press and hold to unlock the wheelchair.



EMERGENCY RELEASE


Locking mechanism can be released manually with the emergency handle on the front edge of the locking plate.



MANUAL RELEASE

Wheelchair can be released manually in the event of vehicle power failure.

ACCESSORIES

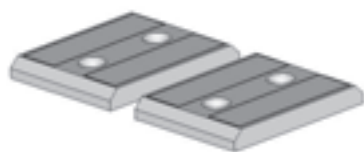

INSTALLATION SPACER

Used to install Permolock on a sloping surface.



The Power of Mobility

Permobil Ltd.
Oceanic, Ground Floor
Navigation Business Park
Waters Meeting Road
Bolton, BL1 8SW
info@permobil.co.uk

www.permobil.com