

permobil

ROHO® HYBRID ELITE
SR™ CUSHION



**THE STABILITY OF FOAM.
THE PROTECTION OF AIR.**

ROHO

THE ROHO HYBRID ELITE SR CUSHION IS DESIGNED FOR PEOPLE WHO:

- have a history of, are at risk for, or currently have any stage pressure injury
- want assistance with setting up their cushion
- have normal or impaired sensation
- require increased sitting stability
- require a more stable environment for transfers
- require lower extremity alignment



FEATURES

Confidence & Independence Available with Smart Check

Take the guesswork out of checking the inflation level of your cushion. Even as your body changes, the optional Smart Check allows you to verify the inflation level to ensure optimal performance and protection.



Positioning

The contoured foam base provides the opportunity for both lower extremity and pelvic alignment and control.



Skin/Soft Tissue Protection

The Hybrid Elite SR cushion hydrostatically loads the ischial tuberosities and the posterior shelf of the greater trochanters for those who are at risk of pressure injuries. The adjustable cushion reduces friction and shear, helping protect you from skin/soft tissue breakdown. It allows for up to 2.75 inches (7 cm) of immersion.



Stability

The contoured foam base provides a solid support under the pelvis and lower extremities to enhance sitting balance and posture. This may also allow for easier transfers in and out of your wheelchair.



SPECIFICATIONS

In the Box, Optional Accessories


Hybrid Elite SR Cushion, two-way stretch cover, foam insert, hand inflation pump, repair kit, and operating instructions. Optional Heavy-Duty cover is available. Smart Check device sold separately.

Construction

ROHO DRY FLOATATION air-filled, adjustable, cellular design insert (cushion insert) made of neoprene (not made with natural rubber latex) with quick disconnect and nickel-plated brass inflation valve. Contoured foam base made of fluid-resistant polyethylene closed-cell foam. Foam insert made of HR polyurethane foam, encased in a water-resistant cover.

Size, Weight, Warranty

Wide range, to fit chair widths 14" to 24"/36 cm to 61 cm
Varies by size, approximately 3 lbs/1.4 kg
(based on cushion size 17" x 18"/43 cm X 45.5 cm)
Cushion: 24-Month Limited Warranty | Smart Check Device: 24-month limited warranty | WEIGHT LIMIT: 500 lbs/225 kg U.S.
MEDICARE CODES: E2622/E2623

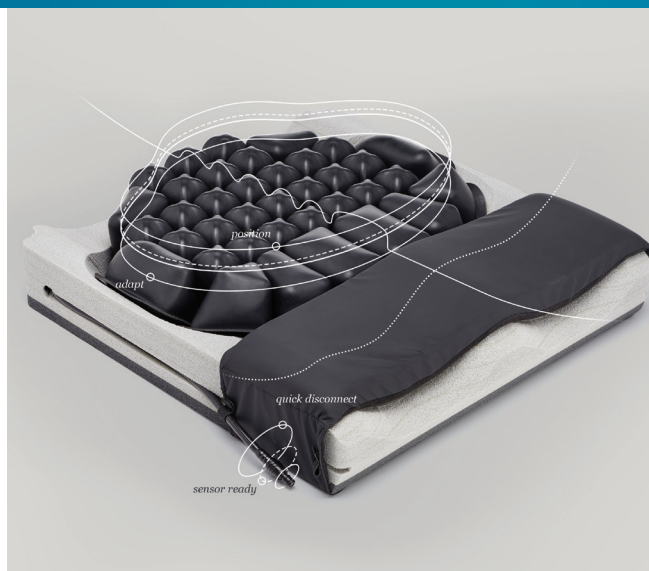

ROHO, Inc.: DRY FLOATATION®, Hybrid Elite SR™, Sensor Ready®, Shape Fitting Technology®, and Smart Check®, and ROHO®.

permobil

ROHO

Permobil Ltd | Office 312, Britannia Office Quarter Waters Meeting, Britannia Way, Bolton, BL2 2HH, England
Tel. (+44) (0)1484 722 888 | info@permobil.co.uk | permobil.co.uk

ROHO, Inc. has a policy of continual product improvement and reserves the right to amend in this document.



The Quality Management System of ROHO, Inc. is certified to ISO 9001:2008 and ISO 13485:2003.

